

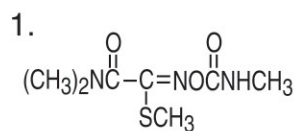


カルバメート系殺虫剤  
Carbamate insecticides


J010220D

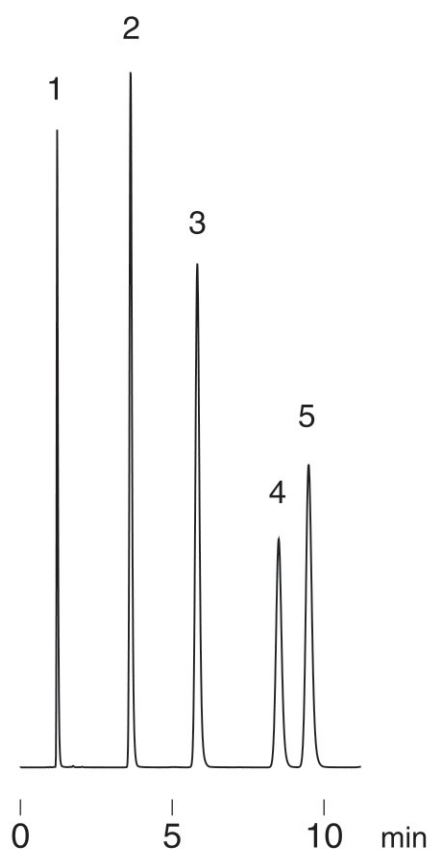


Oxamyl

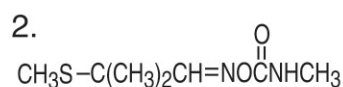


Aldicarb

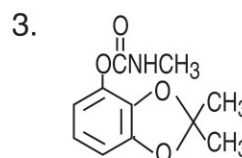

Bendiocarb

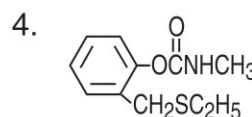

Ethiofencarb

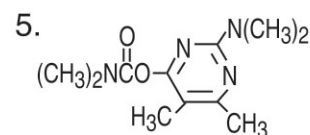

Pirimicarb

Column	: Hydrosphere C18 (3 μm, 12nm) 75 × 4.6 mm I.D.
Eluent	: methanol / 20mM KH <sub>2</sub> PO <sub>4</sub> (45/55)
Flow rate	: 1.0 mL/min
Temperature	: 30 °C
Detection	: UV at 240 nm, 0.64 AUFS
Injection	: 5 μL (0.1~1.5 mg/mL)