


1. Oxytocin (MW 1,007)
2. Met-Enkephalin (MW 574)
3. Neurotensin (MW 1,673)
4. Leu-Enkephalin (MW 556)
5. γ -Endorphin (MW 1,859)
6. β -Endorphin(Human) (MW 3,465)

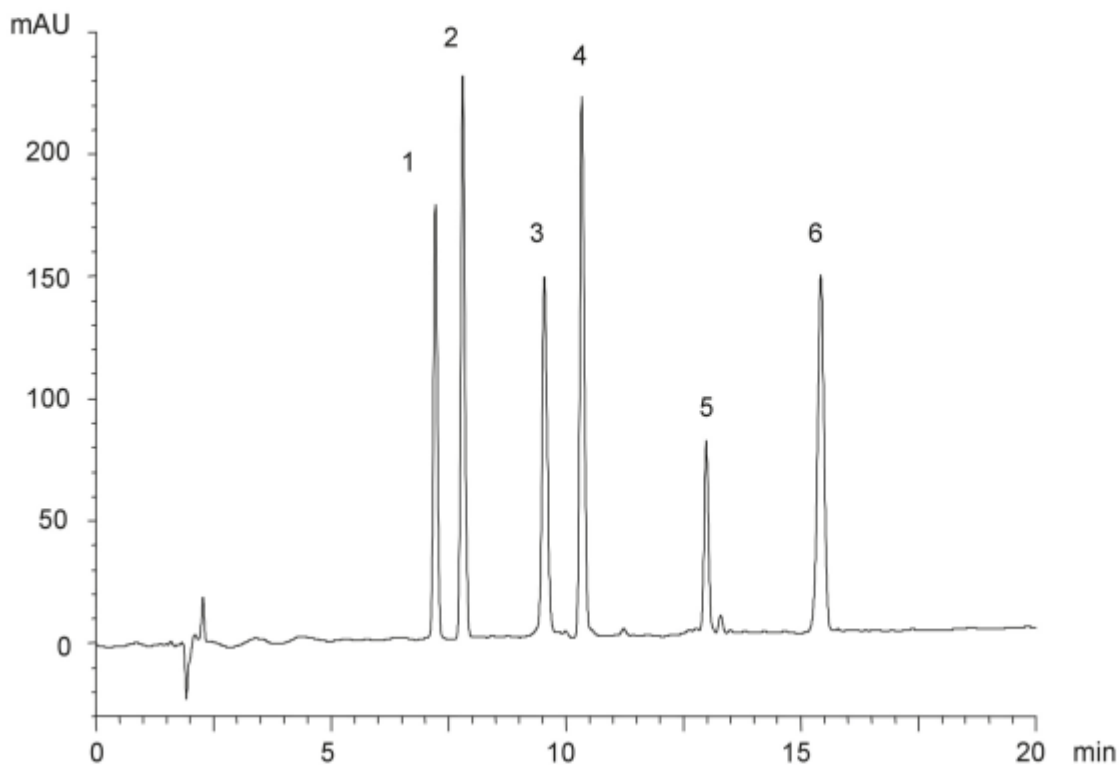

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
 150 \times 4.6 mm I.D.
 Eluent : A) water / TFA (100/0.05)
 B) acetonitrile / TFA (100/0.05)
 20-40%B (0-20min)
 Flow rate : 1.0 mL/min
 Temperature : 37 $^{\circ}$ C
 Detection : UV at 220 nm
 Injection : 10 μ L (0.075~0.25 mg/mL)