

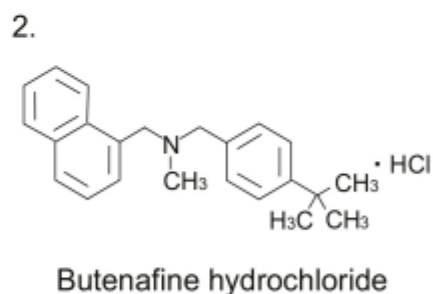
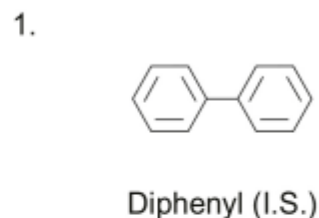


Commercially available reagent
(0.5, 0.17 mg/mL, correspond to standard solution
for butenafine hydrochloride cream)

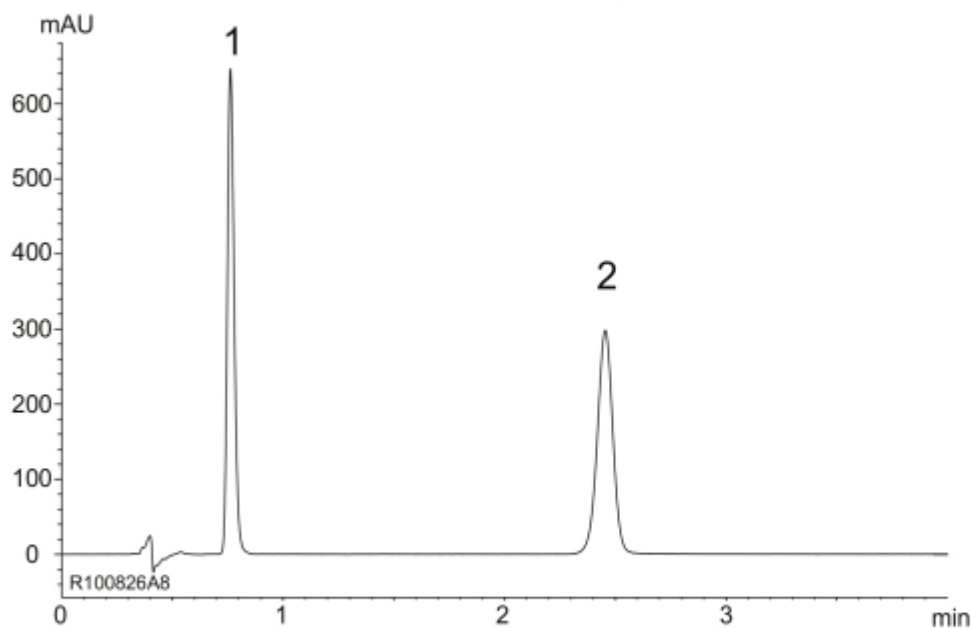

Sample solution for butenafine hydrochloride cream

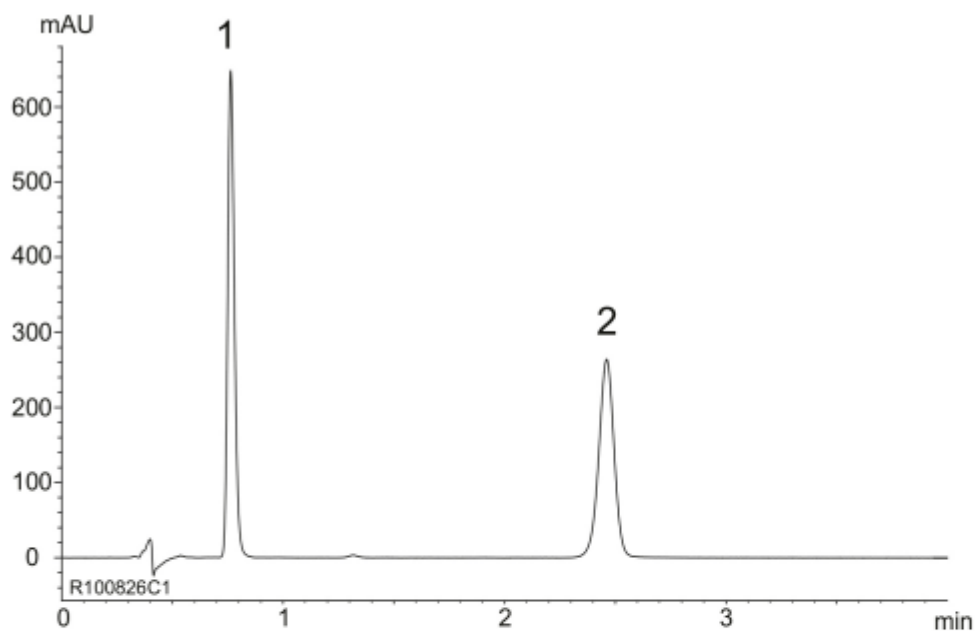

Column : YMC-Pack Pro C8 (3 μ m, 12 nm)
50 X 3.0 mm I.D.
Eluent : 1 mM CH₃COONH₄/acetonitrile (1/4)
Flow rate : 0.7 mL/min (adjust the flow rate so that the retention time of butenafine is about 2.5 min)
Temperature : 40°C
Detection : UV at 282 nm
Injection : 5 μ L

(日本薬局方収載原案記載条件)