
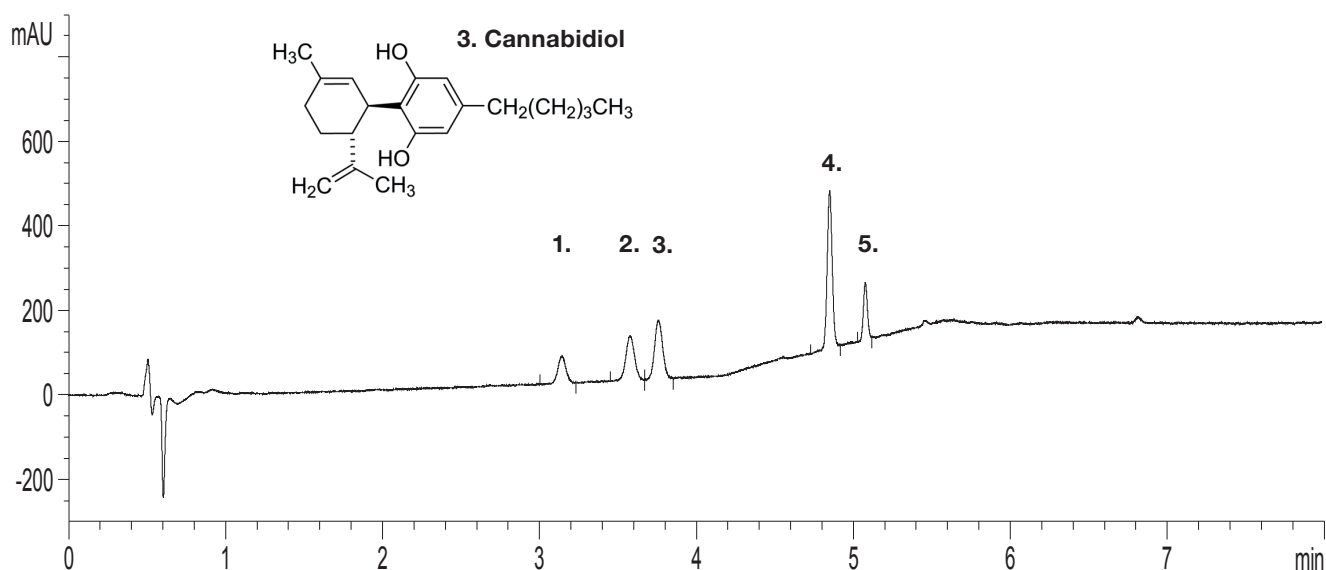
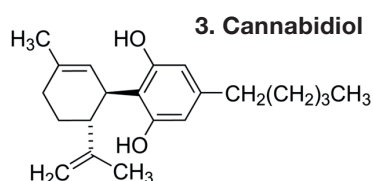
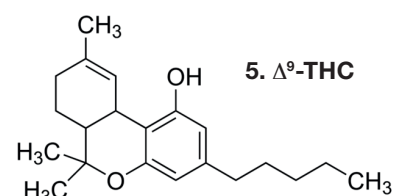
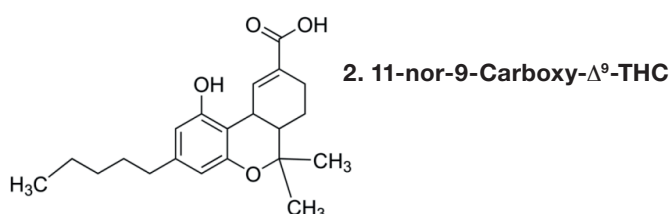
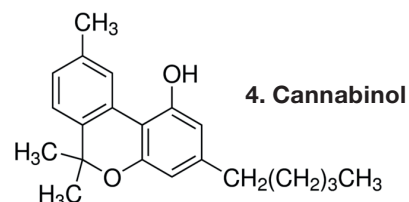


Application Note

Cannabinoids

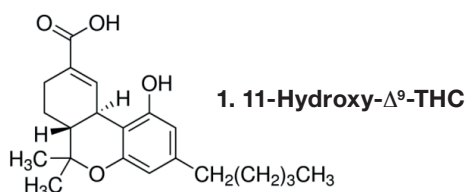

Column: YMC-Triart C18 (1.9 μ m, 12 nm) 100 x 2.0 mm ID
Part No.: TA12SP9-1002PT
Eluent: A) water/formic acid (100/0.1)
B) methanol/formic acid (100/0.1)
Gradient: 75-82 % B (0-3.5 min), 82-100 % B (3.6-5.0 min), 100 % B (5.0-8.0 min)
Flow rate: 0.5 mL/min
Temperature: 25 °C
Detection: UV at 220 nm
Injection: 2 μ L (17 μ g/mL)