



*m*-Cresol

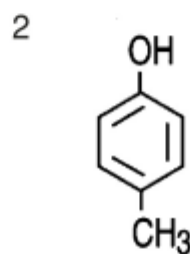

*p*-Cresol

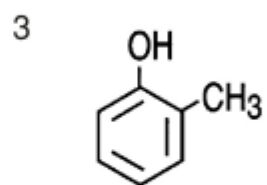

*o*-Cresol

|             |   |
|-------------|---|
| Column      | : J'sphere ODS-H80 (4 $\mu$ m, 8nm)<br>150 $\times$ 4.6 mm I.D. |
| Eluent      | : acetonitrile/water (18/82)                                    |
| Flow rate   | : 1.0 mL/min  |
| Temperature | : 37 $^{\circ}$ C   |
| Detection   | : UV at 260 nm, 0.064 AUFS                                      |
| Injection   | : 10 $\mu$ L (0.3 mg/mL)  |